BAUSCH + LOMB

Bio

true

ONEday (nesofilcon A)

Contact Lenses

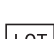
LENS SPECIFIC GUIDE

One contact lens is 22% polymer (nesofilcon A), 78% water and is packaged in a storage solution of sterile borate buffered saline.

 Manufactured by: Bausch & Lomb Incorporated
1400 North Goodman Street, Rochester NY 14609, USA

TM/® are trademarks of Bausch & Lomb Incorporated or its affiliates. © Bausch & Lomb Incorporated.

The following symbols appear on the primary and/or secondary package label:

	Do Not Reuse
	Temperature Limitation
	Sterile Using Steam or Dry Heat
	See Instruction Leaflet
	CE Conformity Marking and the Notified Body Number
	Authorized Representative in European Community
	Federal law restricts this device to sale by or on the order of a licensed practitioner
	Fee Paid for Waste Management
	Use by Date (Expiration Date)
	Batch Code
	Diameter
	Power
	Base Curve

Indication

BAUSCH + LOMB Biotrue ONEday (nesofilcon A) Contact Lens is indicated for the daily wear correction of refractive ametropia (myopia, hyperopia, astigmatism) and reduction of spherical aberration in aphakic and/or non-aphakic persons with non-diseased eyes that exhibit refractive astigmatism up to 2.00 diopters or less. The lens is to be prescribed in spherical powers ranging from +20.000 to -20.000.

Wearing Schedule

The lens is to be prescribed for single-use disposable wear, and is to be discarded after each removal.

Lens Care Procedures

BAUSCH + LOMB Biotrue ONEday (nesofilcon A) Contact Lenses are designed to be disposed of after each daily period of wear. No lens care products are therefore indicated for disinfection of this lens.

Technical Information

Labeled parametric values are determined by specific BAUSCH + LOMB methodology

Manufacture:	Cast Molding
Material:	Nesofilcon A (Non-Ionic)
Water Content:	78%
Dk:	42
Sphere Powers:	+6.50 to -6.50D (in 0.25D steps) -7.00 to -9.00D (in 0.50D steps)
Diameter:	14.2mm
Centre Thickness:	0.080mm to 0.240mm
Base Curve:	8.6mm
Recommended Usage Modality:	Daily Disposable

RECOMMENDED TRIAL FITTING PROCEDURE FOR Biotrue ONEday (nesofilcon A) Contact Lenses

Selection of a Trial Lens

- Carry out the preliminary examination in order to assess the patient's suitability for contact lens wear
- Select the initial Biotrue ONEday trial lens as follows:
- The initial lens power is predicted by the spherical equivalent refraction.

Trial Fitting

- Allow the lens to settle on the eye for 10-20 minutes. The time required will vary with the individual.
- Wait for any increase in tear flow to subside before evaluating the lens fit.

Evaluation

- If patients demonstrate reduced or unstable visual acuity, excessive movement in primary gaze or up gaze, poor corneal coverage, limbal indentation or any signs of a poor fitting lens, remove the trial lens and try an alternative design.

Dispensing

Determine Final Lens Power :

Dispense :

- Lens
- Instructions For Use

Discuss :

- Wearing schedule
- Follow-up schedule
- Adaptive symptoms
- Warning signs

GENERAL GUIDE FOR BAUSCH + LOMB CONTACT LENSES

Important: This information should be explained to the patient when supplied with BAUSCH + LOMB Contact Lenses

Intended Use

BAUSCH + LOMB Contact Lenses can be prescribed for the correction of vision defects including myopia, hyperopia, astigmatism, presbyopia, in aphakia and/or non-aphakia persons.

Lens Care Procedures

BAUSCH + LOMB recommends that BAUSCH + LOMB lens care products be used with BAUSCH + LOMB Contact Lenses. Please refer to the appropriate guide for details. The practitioner should recommend a care system to the patient for these lenses. For whatever system is selected, the instructions which accompany these care products must be explained to the patient.

Contraindications

Use your professional judgement and experience in evaluating patient or ocular suitability for contact lenses. In addition, BAUSCH + LOMB lenses are contraindicated when the following conditions exist:

- Unhealthy patient eyes
- History of patient non-compliance with contact lens care and disinfection regimens, wearing restrictions, wearing schedule, or follow-up visit schedule.
- Patient inability or unwillingness, because of age, infirmity or other mental or physical conditions, or an adverse working or living environment, to understand or comply with any warnings, precautions, restrictions, or directions.

Warnings

After a thorough eye examination, including appropriate medical background, patients should be fully apprised by the prescribing practitioner of all the risks with contact lens wear. Patients should be advised of the following warnings pertaining to contact lens wear:

- Problems with contact lenses and lens care products could result in **serious injury** to the eye. It is essential that patients follow their eye care practitioner's direction and all labeling instructions for proper use of lenses and lens care products, including the lens case. Eye problems, including corneal ulcers, can develop rapidly and lead to **loss of vision**.
- When prescribed for the Frequent Replacement Program, the need for strict compliance with the care regimen including cleaning of the lens case, wearing restrictions, wearing schedule and follow-up visit schedule should be emphasized to the patient.
- Studies have shown that contact lens wearers who are smokers have a higher incidence of adverse reactions than non-smokers.

Extended Wear

- The risk of microbial keratitis has been shown to be greater among users of extended wear soft hydrophilic contact lenses than among users of daily wear contact lenses. The risk among these extended wear lens users increases with the number of consecutive days that the lenses are worn between removals, beginning with the first overnight use. Some researchers believe that these complications are caused by one or more of the following: a weakening of the cornea's resistance to infections, particularly during a closed-eye condition, as a result of hypoxia; an eye environment which is somewhat more conducive to the growth of bacteria and other microorganisms, particularly when a regular periodic lens removal and disinfecting or disposal schedule has not been adhered to by the patient; improper lens disinfection or cleaning by the patient; contamination of lens care products; poor personal hygiene by the patient; patient unsuitability to the particular lens or wearing schedule; accumulation of lens deposits; damage to the lens; improper fitting; length of wearing time; and the presence of ocular debris or environmental contaminants. While the great majority of patients successfully wear contact lenses, extended wear of lenses also is reported to be associated with a higher incidence and degree of epithelial microcysts and infiltrates, and endothelial polymegathism, which require consideration of discontinuation or restriction of extended wear. The epithelial conditions are reversible upon discontinuation of extended wear.

The reversibility of endothelial effects of contact lens wear has not been conclusively established. As a result, practitioners views of extended wearing times have varied until now from not prescribing extended wear at all to prescribing flexible wearing times from occasional overnight wear to prescribing extended wearing periods from 1 to 7 days with specified intervals of no lens wear for certain patients, with follow-up visits, and with proper care regimen.

New silicone hydrogel lenses such as PureVision®, PureVision® Toric, PureVision® Multi-Focal and PureVision2® offer benefits over conventional hydrogel lenses and may be fitted for up to 30 days of continuous wear. Clinical data to date suggests that these lenses provide fewer adverse effects than conventional hydrogel lenses for extended wear but a post-marketing surveillance study has not been completed, and no conclusions can be drawn about whether the risk of microbial keratitis in extended wear is reduced with silicone hydrogel lenses.

- If a patient experiences eye discomfort, excessive tearing, vision changes, or redness of the eye, the patient should be instructed to **immediately remove lenses** and promptly contact his or her eye care practitioner.

Additional Warnings and Precautions for Biotrue ONEday (nesofilcon A) Contact Lenses:

Warning: UV absorbing contact lenses are NOT substitutes for protective UV absorbing eyewear such as UV absorbing goggles or sunglasses because they do not completely cover the eye and surrounding area. You should continue to use UV absorbing eyewear as directed.

Note: Long term exposure to UV radiation is one of the risk factors associated with cataracts. Exposure is based on a number of factors such as environmental conditions (altitude, geography, cloud cover) and personal factors (extent and nature of outdoor activities). UV blocking contact lenses help provide protection against harmful UV radiation.

Note: The effectiveness of wearing UV-absorbing contact lenses in preventing or reducing the incidence of ocular disorders associated with exposure to UV-light has not been established at this time. However, clinical studies have not been done to demonstrate that wearing UV blocking contact lenses reduce the risk of developing cataracts or other eye disorders. Consult your Eye Care Professional for more information.

Precautions

- Contact lens wear may not be suitable for certain occupations, or, in other instances, may require eye protection equipment.
- Environmental fumes, smoke, dust, vapors and windy conditions must be avoided in order to minimize the chances of lens contamination or physical trauma to the cornea.
- Hard contact lens solutions, not indicated for use with soft lenses, may not be used in the soft lens care system. Soft lens solutions should not be used with hard lenses. Serious corneal injury may result from mismatching lens and solutions types.
- Eye injury from irritation or infection and damage to lenses may result if cosmetics, lotions, soaps, creams, hair sprays or deodorants come into contact with lenses.
- Tweezers or other tools should not be used by patients to remove a lens from the lens container. The lens should be poured into the hand or removed with a finger.
- Patients must be instructed on and demonstrate the ability to promptly remove the lenses from their eyes.
- Fluorescein should not be used while SOFT lenses are in the patient's eye. Soft lenses absorb this dye and become discolored. Fluorescein in the eye should be thoroughly flushed with a sterile saline solution and the soft lens should not be reinserted until at least one hour has elapsed.

Adverse Effects

The following symptoms may occur:

- Eye pain
- Eye stinging, burn or itch (irritation)
- Comfort is less than when lens was first placed on eye
- Feeling of something in the eye (foreign body, scratched area)
- Excessive watering (tearing) of the eye
- Unusual eye secretions
- Redness of the eyes
- Reduced sharpness of vision (poor visual acuity)
- Blurred vision, rainbows or haloes around objects
- Change in sensitivity to light (photophobia)
- Feeling of dryness

The patient should be instructed that if any of the above symptoms occur they should:

- Immediately remove the lenses
- If the discomfort or problems stop, then look closely at the lens
- If the lens is in any way damaged, do not put the lens back in the eye. Place the lens in the storage case and contact the eye care practitioner.
- If the lens has dirt, an eyelash or other foreign body on it, or the problem stops and the lens appears undamaged, thoroughly clean, rinse and disinfect the lenses, then reinsert them.
- If the above symptoms continue after removal of the lens or recur upon reinsertion of the lens, the lenses should be removed immediately and the patient should immediately contact the eye care practitioner or a physician, who must determine the need for examination, treatment or referral without delay.

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